

Primer Design Ltd

Z-Path-HRV16

Human Rhinovirus 16

Kit version: 1

Target region:

Polyprotein gene

genesig[®] Advanced Kit

150 tests

GENESIG

Kits by Primerdesign

For general laboratory and research use only

Introduction to Human Rhinovirus 16

Rhinovirus is a genus of the Picornaviridae family of viruses. Rhinoviruses are the most common viral infective agents in humans and are a causative agent of the common cold, being responsible for approximately 50% of all cases. Rhinoviruses have linear, single stranded, positive-sense RNA genomes of between 7.2 and 8.5 kb in length. At the 5' end of the genome is a virus-encoded protein, and like mammalian mRNA, there is a 3' poly-A tail. Structural proteins are encoded in the 5' region of the genome and non structural at the 3' end. The viral particles themselves are not enveloped and are icosahedral in structure with a diameter of approximately 30nm.

Rhinoviruses are composed of a capsid that contains four viral proteins VP1, VP2, VP3 and VP4. VP1, VP2, and VP3 form the major part of the protein capsid. The much smaller VP4 protein has a more extended structure and lies at interface between the capsid and the RNA genome. There are 60 copies of each of these proteins assembled as an icosahedron. Antibodies are a major defence against infection with the epitopes lying on the exterior regions of VP1-VP3.

Rhinoviruses have two main modes of transmission: direct transmission via inhalation of aerosolised respiratory droplets and indirect transmission via contact with surfaces contaminated with respiratory droplets. The viral capsid proteins bind to specific nasal cell receptors mediating entry into the host cell. Intracellular enzymes then digest the viral capsid leaving only the genome which functions as a template for the synthesis of the complementary strand. This encodes the viral capsid proteins for new viral particles as well as being packaged into these. The newly synthesised virus buds off from the host cell membrane facilitating infection of other cells.

Rhinoviruses occur worldwide causing disease with symptoms such as fever, cough, and nasal congestion which are the result of the immune system's response to the virus. The frequency of colds is high in childhood and decreases during adulthood most probably because of the possession of immunity.

There are many HRV serotypes which, based on differences in their surface proteins and genome sequences have been segregated into two different groups: Human Rhinovirus A and B although some serotypes remain unclassified. HRV16 belongs to the HRVA species.

Specificity

The genesig® Advanced Kit for Human Rhinovirus 16 (HRV16) is designed for the in vitro detection and quantification of HRV16 genomes. The kit is designed to have a broad detection profile. Specifically, the primers will detect over 95% of sequences available on the NCBI database at the time of last review.

The dynamics of genetic variation mean that new sequence information may become available after the most recent review. If you require further information or have a specific question about the detection profile of this kit then please send an e-mail to techsupport@primerdesign.co.uk and our team will answer your question.

Kit contents

- **1x HRV16 primer/probe mix (150 reactions, BROWN)**
FAM labelled
- **1x HRV16 positive control template (for Standard curve, RED)**
- **1x Internal extraction control primer/probe mix (150 reactions, BROWN)**
VIC labelled as standard
- **1x Internal extraction control RNA (150 reactions, BLUE)**
- **1x Endogenous control primer/probe mix (150 reactions, BROWN)**
FAM labelled, Target: Human ACTB gene as standard
- **1x RNase/DNase free water (WHITE)**
for resuspension of primer/probe mixes
- **3x Template preparation buffer (YELLOW)**
for resuspension of internal control template, positive control template and standard curve preparation

Reagents and equipment to be supplied by the user

Real-time PCR instrument

Extraction kit

This kit is recommended for use with genesig® Easy DNA/RNA Extraction kit or exsig®Mag. However, it is designed to work well with all processes that yield high-quality nucleic acid with minimal PCR inhibitors.

oasig® lyophilised OneStep or Precision Plus® OneStep 2X RT-qPCR Master Mix

This kit is intended for use with oasig® lyophilised OneStep or PrecisionPLUS® OneStep 2X RT-qPCR Master Mix.

Pipettors and filter tips

Vortex and centrifuge

1.5 ml microtubes

qPCR plates or reaction tubes

Kit storage and stability

This kit is stable for shipping at ambient temperature but should be stored at -20°C upon arrival. Once the lyophilised components have been resuspended, they should not be exposed to temperatures above -20°C for longer than 30 minutes at a time and unnecessary repeated freeze/thawing should be avoided. The kit is stable for six months from the date of resuspension under these circumstances.

If a standard curve dilution series is prepared this can be stored frozen for an extended period. If you see any degradation in this serial dilution a fresh standard curve can be prepared from the positive control. Primer Design Ltd does not recommend using the kit after the expiry date stated on the pack.

Suitable sample material

This kit can be used with all types of samples from various origins. Please ensure that the extracted nucleic acid samples are suitable in terms of purity, concentration, and RNA integrity. An internal extraction control is provided to test for non-specific PCR inhibitors.

Dynamic range of test

Under optimal PCR conditions the kit can achieve priming efficiencies between 90-110% and detect less than 100 copies of target template. If running a positive control standard curve for a quantitative result, and an efficiency of between 90% to 110% is not achieved, then the run should be repeated with a freshly prepared standard curve.

Principles of the test

Real-time PCR

A target specific primer and probe mix is provided, and this can be detected through the FAM channel.

The primer and probe mix provided exploits the so-called TaqMan® principle. During PCR amplification, forward and reverse primers hybridize to the target cDNA. A fluorogenic probe is included in the same reaction mixture which consists of a DNA probe labelled with a 5'-dye and a 3'-quencher. During PCR amplification, the probe is cleaved, and the reporter dye and quencher are separated. The resulting increase in fluorescence can be detected on a range of qPCR platforms.

Positive control

For copy number determination and as a positive control for the PCR set up, the kit contains a positive control template. This can be used to generate a standard curve of target copy number/Cq value. Alternatively, the positive control can be used at a single dilution where full quantitative analysis of the samples is not required. Each time the kit is used, at least one positive control reaction must be included in the run. A positive result indicates that the primers and probes for detecting the target gene worked properly in that particular experimental scenario. If a negative result is obtained the test results are invalid and must be repeated. Care should be taken to ensure that the positive control does not contaminate any other kit component which would lead to false-positive results. This can be achieved by handling this component in a Post PCR environment. Care should also be taken to avoid cross-contamination of other samples when adding the positive control to the run. This can be avoided by sealing all other samples and negative controls before pipetting the positive control into the positive control well.

Negative control

To validate any positive findings a negative control reaction should be included every time the kit is used. For this reaction the RNase/DNase free water should be used instead of template. A negative result indicates that the reagents have not become contaminated while setting up the run. It is also known as a No Template Control or NTC.

Internal RNA extraction control

When performing RNA extraction, it is often advantageous to have an exogenous source of RNA template that is spiked into the lysis buffer. This control RNA is then co-purified with the sample RNA and can be detected as a positive control for the extraction process. Successful co-purification and qPCR for the control RNA also indicates that PCR inhibitors are not present at a high concentration.

A separate qPCR primer/probe mix is supplied with this kit to detect the exogenous RNA using qPCR. The PCR primers are present at PCR limiting concentrations which allows multiplexing with the target sequence primers. Amplification of the control RNA does not interfere with detection of the target RNA even when present at low copy number. The Internal control is detected through the VIC channel and gives a Cq value of 28+/-3 depending on the level of sample dilution.

Endogenous control

A primer/probe mix for detection of the endogenous control gene is included in the kit, which allows confirmation of a valid biological sample from this host. Detection of the endogenous control is through the FAM channel, and it is therefore NOT possible to perform a multiplex reaction with the target specific primer/probe mix. Amplification of the endogenous control may depend on the sample type used. Please note that if samples from a different species are used, the endogenous control may not be appropriate, but the internal extraction control is advised to be used.

Resuspension protocol

To minimize the risk of contamination with foreign RNA/DNA, we recommend that all pipetting is performed in a PCR clean environment. Ideally this would be a designated PCR lab or PCR cabinet. Filter tips are recommended for all pipetting steps.

1. Pulse-spin each tube in a centrifuge before opening.

This will ensure lyophilised primer/probe mix or template is in the base of the tube and is not lost upon opening the tube.

2. Resuspend the primer/probe mixes in the RNase/DNase free water supplied, according to the table below:

To ensure complete resuspension, vortex each tube thoroughly.

Component - resuspend in water	Volume
Pre-PCR pack	
HRV16 primer/probe mix (BROWN)	165 µl
Internal extraction control primer/probe mix (BROWN)	165 µl
Endogenous control primer/probe mix (BROWN)	165 µl

3. Resuspend the internal control template and positive control template in the template preparation buffer supplied, according to the table below:

To ensure complete resuspension, vortex each tube thoroughly.

Component - resuspend in template preparation buffer	Volume
Pre-PCR heat-sealed foil	
Internal extraction control RNA (BLUE)	600 µl
Post-PCR heat-sealed foil	
HRV16 Positive Control Template (RED) *	500 µl

* This component contains high copy number template and is a VERY significant contamination risk. It must be opened and handled in a separate laboratory environment, away from the other components.

RNA extraction

The internal extraction control RNA can be added either to the RNA lysis/extraction buffer or to the RNA sample once it has been resuspended in lysis buffer.

DO NOT add the internal extraction control RNA directly to the unprocessed biological sample as this will lead to degradation and a loss in signal.

- 1. Add 4µl of the Internal extraction control RNA (**BLUE**) to each sample in RNA lysis/extraction buffer.**
- 2. Complete RNA extraction according to the manufacturer's protocols.**

OneStep RT-qPCR detection protocol

For optimum performance and sensitivity.

All pipetting steps and experimental plate set up should be performed on ice. After the plate is prepared proceed immediately to the OneStep amplification protocol. Prolonged incubation of reaction mixes at room temperature can lead to PCR artifacts that reduce the sensitivity of detection.

1. For each RNA sample prepare a reaction mix according to the table below: Include sufficient reactions for positive and negative controls.

Component	Volume
oasig [®] OneStep or Precision [®] PLUS OneStep 2X RT-qPCR Master Mix	10 µl
HRV16 primer/probe mix (BROWN)	1 µl
Internal extraction control primer/probe mix (BROWN)	1 µl
RNase/DNase free water (WHITE)	3 µl
Final Volume	15 µl

2. (Optional) For each RNA sample prepare an endogenous control reaction according to the table below.

This control reaction will provide useful information regarding the quality of the biological sample.

Component	Volume
oasig [®] OneStep or Precision [®] PLUS OneStep 2X RT-qPCR Master Mix	10 µl
Endogenous control primer/probe mix (BROWN)	1 µl
RNase/DNase free water (WHITE)	4 µl
Final Volume	15 µl

3. Pipette 15µl of these mixes into each well according to your experimental qPCR plate set-up.
4. Pipette 5µl of RNA template into each well according to your experimental plate set-up.

For negative control wells use 5µl of RNase/DNase free water (**WHITE**). For positive control wells use 5µl of the positive control template (**RED**). The final volume in each well is 20µl.

5. (Optional) Standard curve preparation for quantitative analysis.

For quantitative analysis of the samples a standard curve dilution series can be prepared using the positive control template (**RED**). This is not required for qualitative analysis.

5.1 Reaction mix preparation for the standard curve.

Include sufficient reactions for each dilution of the standard curve:

Component	Volume
oasig® OneStep or Precision®PLUS OneStep 2X RT-qPCR Master Mix	10 µl
HRV16 primer/probe mix (BROWN)	1 µl
RNase/DNase free water (WHITE)	4 µl
Final Volume	15 µl

5.2 Preparation of a 10-fold standard curve dilution series.

- Pipette 90 µl of template preparation buffer (**YELLOW**) into 5 tubes and label them tube 2-6. The neat positive control tube (**RED**) is considered tube 1.
- Pipette 10 µl of positive control template (**RED**) into tube 2.
- Vortex thoroughly.
- Change pipette tip and pipette 10 µl from tube 2 into tube 3.
- Vortex thoroughly

Repeat steps **d** and **e** across the tubes to complete the dilution series.

Standard Curve	Copy Number
Tube 1 Positive control (RED)	2 x 10 ⁵ per µl
Tube 2	2 x 10 ⁴ per µl
Tube 3	2 x 10 ³ per µl
Tube 4	2 x 10 ² per µl
Tube 5	20 per µl
Tube 6	2 per µl

5.3 Pipette 15µl of reaction mix and 5µl of the respective standard into each well for the standard curve according to your plate set-up.

The final volume in each well is 20µl.

OneStep RT-qPCR Amplification Protocol

Recommended amplification conditions when using oasisig® OneStep or Precision®PLUS OneStep 2X RT-qPCR Master Mix.


	Step	Time	Temp
	Reverse Transcription	10 min	55 °C
	Enzyme activation	2 min	95 °C
Cycling x50	Denaturation	10 s	95 °C
	DATA COLLECTION *	60 s	60 °C

* Fluorogenic data should be collected during this step through the FAM and VIC channels


Interpretation of results

Target (FAM)	Internal control (VIC)	Positive control	Negative control	Interpretation
≤ 30	+ / -	+	-	POSITIVE QUANTITATIVE RESULT calculate copy number
> 30	+	+	-	POSITIVE QUANTITATIVE RESULT calculate copy number
> 30	-	+	-	POSITIVE QUALITATIVE RESULT do not report copy number as this may be due to poor sample extraction
-	+	+	-	NEGATIVE RESULT
+ / -	+ / -	+	≤ 35	EXPERIMENT FAILED due to test contamination
+ / -	+ / -	+	> 35	*
-	-	+	-	SAMPLE PREPARATION FAILED
+ / -	+ / -	-	+ / -	EXPERIMENT FAILED

*Where the test sample is positive and the negative control is positive with a Cq > 35 , the sample must be reinterpreted based on the relative signal strength of the two results:

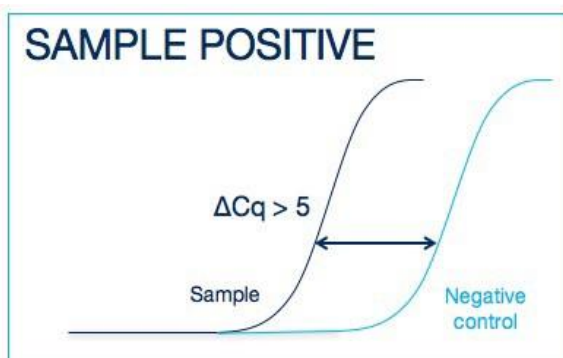

If the sample amplifies > 5 Cq earlier than the negative control, then the sample should be reinterpreted (via the table above) with the negative control verified as negative.

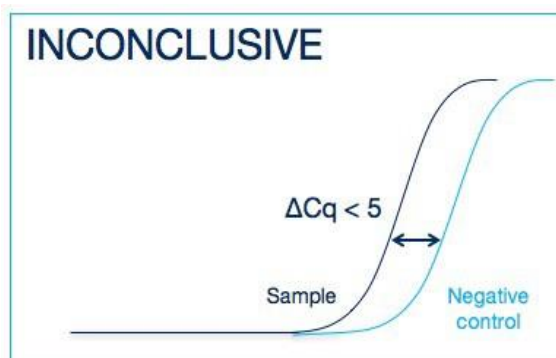

If the sample amplifies < 5 Cq earlier than the negative control, then the positive sample result is invalidated, and the result should be determined inconclusive due to test contamination. The test for this sample should be repeated.

Positive Control

The positive control template is expected to amplify between Cq 16 – 23 in the FAM channel. Failure to satisfy this quality control criterion is a strong indication that the experiment has been compromised and should be repeated.

Internal PCR control

The Cq value obtained with the internal control will vary significantly depending on the extraction efficiency, the quantity of RNA added to the RT and PCR reaction and the individual machine settings. Cq values of 28±3 are within the normal range. When amplifying a sample with a high genome copy number, the internal extraction control may not produce an amplification plot. This does not invalidate the test and should be interpreted as a positive experimental result.

Endogenous control

The signal obtained from the endogenous control primer and probe set will vary according to the amount of biological material present in a given sample. An early signal indicates the presence of a good yield of biological material. A late signal suggests that little biological material is present in the sample.

Please note that if samples from a different species or environmental samples are used, the endogenous control may not be appropriate and would not result in amplification.

Notices and disclaimers

This product is developed, designed and sold for research purposes only. It is not intended for human diagnostic or drug purposes or to be administered to humans unless clearly expressed for that purpose by the Food and Drug Administration in the USA or the appropriate regulatory authorities in the country of use. During the warranty period Primer Design Ltd genesig® detection kits allow precise and reproducible data recovery combined with excellent sensitivity. For data obtained by violation to the general GLP guidelines and the manufacturer's recommendations the right to claim under guarantee is expired. PCR is a proprietary technology covered by several US and foreign patents. These patents are owned by Roche Molecular Systems Inc. and have been sub-licensed by PE Corporation in certain fields. Depending on your specific application you may need a license from Roche or PE to practice PCR. Additional information on purchasing licenses to practice the PCR process may be obtained by contacting the Director of Licensing at Roche Molecular Systems, 1145 Atlantic Avenue, Alameda, CA 94501 or Applied Biosystems business group of the Applied Biosystems Corporation, 850 Lincoln Centre Drive, Foster City, CA 94404. In addition, the 5' nuclease assay and other homogeneous amplification methods used in connection with the PCR process may be covered by U.S. Patents 5,210,015 and 5,487,972, owned by Roche Molecular Systems, Inc, and by U.S. Patent 5,538,848, owned by The Perkin-Elmer Corporation.

Trademarks

PrecisionPLUS® is a trademark of Primer Design Ltd.

genesig® is a registered trademark of Primer Design Ltd.

oasig® is a registered trademark of Primer Design Ltd.

exsig®Mag is a registered trademark of Primer Design Ltd.

The PCR process is covered by US Patents 4,683,195, and 4,683,202 and foreign equivalents owned by Hoffmann-La Roche AG. TaqMan® is a registered trademark of Roche Molecular Systems, Inc., The purchase of the Primer Design Ltd reagents cannot be construed as an authorization or implicit license to practice PCR under any patents held by Hoffmann-La Roche Inc.